



# VETROSON® MILLENNIUM 25/30 ULTRASONIC SCALER

## STATE OF THE ART - MAKES OTHER SCALERS OBSOLETE

The VETROSON® Millennium 25/30 Ultrasonic Scaler accepts the new, superior 30k inserts as well as the standard 25k inserts.

VETROSON® Millennium 25/30 Ultrasonic Scaler features a new concept in autotuning. The unit zeros in on an insert's natural resonance, producing a better and more consistent spray pattern for a broader spectrum of inserts.

The VETROSON® Millennium 25/30 Ultrasonic Scaler may be adjusted to a "low" power setting for perio work or "full" power for scaling above the gums with excellent cavitation and spraying patterns at both settings. Current scalers with one power setting cannot be used for working under the gums.

VETROSON® Millennium 25/30 Ultrasonic Scaler has a power booster system. This automatically increases the power to a point midway between the perio setting and maximum power. Push the button again to resume normal subgingival, debridement.

VETROSON® Millennium 25/30 Ultrasonic Scaler has a detachable, heat resistant handpiece sheath which twists off for autoclaving. This feature is patented.

***New Perio Insert and Perio Setting allow work under gums without damaging tissue.***

A separate transformer generates a minimum of heat.

**THREE YEAR WARRANTY**

**Used and Recommended by Most Veterinary Board Certified Dentists**



The basic unit includes a VETROSON® "SS" #10 Insert 30k, a VETROSON "SS" Universal Perio Insert 30k, Filter & Quick Disconnect water hook up.



It's compact!  
1 3/4" x 5 1/2" x 7" weighing 2 3/4 lbs.



**SUMMIT HILL LABORATORIES**  
Tinton Falls Business Center  
One Sheila Drive  
Tinton Falls, New Jersey 07724

# VETROSON® “SS” (STAINLESS STEEL) INSERTS

**DISTINCTLY SUPERIOR**


## FOR VETROSON® & CAVITRON® TYPE ULTRASONIC SCALERS

The new inserts are made of stainless steel as opposed to traditional scaler inserts which are a combination of plastic and metal. The result of stainless steel uni-construction is a more durable scaling instrument, which will not crack or loosen upon repeated autoclaving.




## THE NEW VETROSON® STAINLESS STEEL “SS” INSERT

Inserts with the usual #10 tip configurations, are available in both 25k (for present machines) and 30k (for the new **VETROSON® MILLENNIUM 25/30 ULTRASONIC SCALER**). The new 40% thinner universal Perio insert is only available in 30k. A color coded “O” ring will identify the inserts’ operating frequencies, black for 25k and green for 30k inserts.

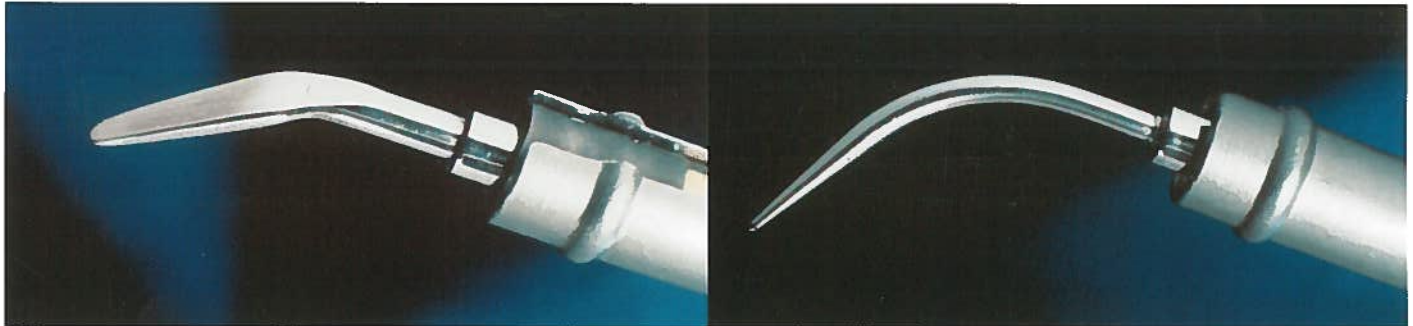

**VETROSON® “SS” #3 INSERT**

**VETROSON® “SS” #10 INSERT**

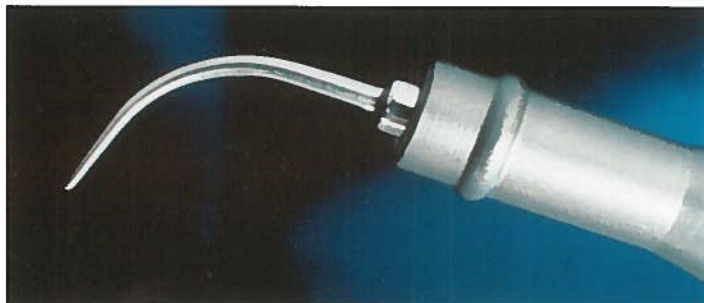

**VETROSON® “SS” PERIO INSERT - 40% THINNER**